BLUEPRINT proposed data submission schemas Draft v0.2.3.2-566596e769717c140078683abfe46f007f199452

(model SHA1 a4667cf8d285d0b4339bd510c4492dfcf2775b7a)

(external controlled vocabulary SHA1 5fc45065b5d8fbe1e94c3370acf8bdd965540acc)

José María Fernández

March 13, 2013

Contents

1	Data	Submission	1
	1.1	Overview of Data Submission Process	1
	1.2	Preparing Sample Tracking Data and Analyzed Contents for their submission	1
		1.2.1 File Naming Conventions	2
		1.2.2 Tabular File Structure	3
	1.3	File Submission Procedure	4
2	DCC	Submission Tabular Formats	5
_	2.1	Sample Tracking Submission File Specifications	5
	2.1	2.1.1 Donor Data File	6
		2.1.2 Specimen Data File	6
		2.1.3 Analyzed Sample Data File	8
	2.2	Gene Expression	9
	2.2	2.2.1 Expression - Metadata File	9
		2.2.2 Expression - Gene File	10
	2.3		12
	2.5		12
			15
	2.4		17
	۷.٦		17
			18
		g g	19
	2.5		20
	2.5		20
			21
			22
	2.6		23
	2.0		23
			24
			26
		2.0.5 DIVIV -tatton - Secondary / matysts i tte	20
Α	Cont	rolled Vocabularies	27
	A.1	Ensembl Genes	27
	A.2	Ensembl Transcripts	27
	A.3		27
	A.4		27
	A.5		27
	A.6	Value Codes for Platform or Validation Platform	28
	A.7	Chromosome Names for Reference Genome GRCh37	30
	8.A	Value Codes for Reference Genome Assembly Version	33
	A.9	Value Codes for Raw Data Repository	33

List of Figures

2.1	Sample Tracking Data Sub-Schema	-
2.2	Gene Expression Sub-Schema	Ć
2.3	Junction Read Count explanation	2
2.4	Exon Junction Sub-Schema	2
2.5	Protein-DNA interactions Sub-Schema	7
2.6	Regulatory Regions Sub-Schema	2(
2.7	DNA Methylation , Hydroxy-Methylation, Formylation, etc Sub-Schema	23



List of Tables

2.1.1 Donor Data File
2.1.2 Specimen Data File
2.1.3 Analyzed Sample Data File
2.2.1 Expression - Metadata File
2.2.2 Expression - Gene File
2.3.1 Exon Junction - Metadata File
2.3.2 Exon Junction - Primary Analysis File
2.4.1 Protein-DNA interaction - Metadata File
2.4.2 Protein-DNA interaction - Primary Analysis File
2.4.3 Protein-DNA interaction - Secondary Analysis File
2.5.1 Regulatory regions - Metadata File
2.5.2 Regulatory regions - Primary Analysis File
2.5.3 Regulatory regions - Secondary Analysis File
2.6.1 DNA *-lation - Metadata File
2.6.2 DNA *-lation - Primary Analysis File
2.6.3 DNA *-lation - Secondary Analysis File
A.4.1 ISO 3166-1 and ISO 3166-2 aliases
A.5.1 Institution ID
A.6.1 Value Codes for Platform or Validation Platform
A.7.1 Chromosome Names for Reference Genome GRCh37
A.8.1 Value Codes for Reference Genome Assembly Version
A.9.1 Value Codes for Raw Data Repository

Chapter 1

Data Submission

1.1 Overview of Data Submission Process

There are four major steps in the data submission process:

- 1. Submit raw sequence data to the European Genome-phenome Archive
- 2. Prepare the BLUEPRINT submission files according to DCC data format specifications
- 3. Verify conformity of the submission files
- 4. Submit files to the DCC Secure FTP server

All submitted data must be based on Human reference genome assembly GRCh37 and Ensembl gene set version 70

When submitting experimental data, please make sure you've already deposited your raw data to the appropriate public data repositories (eg: sequencing reads to EBI EGA) and then populate in your submission files the data elements **raw_data_repository** and **raw_data_accession** with the correct repository and accession number respectively.

1.2 Preparing Sample Tracking Data and Analyzed Contents for their submission

Submitted experimental data files must be from any one of these categories:

- Sample Tracking
- Gene Expression
- Exon Junctions
- DNA *-lation (Methylation, Hydroxy-Methylation, Formylation, etc...)
- Protein-DNA interactions
- Regulatory regions

BLUEPRINT DCC is hosting both sample tracking data and analyzed contents. Contents must be sent following the textual tabular formats defined below. Files with those contents must also follow the BLUEPRINT DCC file naming convention.

Each submitter must have a unique signing key, provided by DACO and DCC. Each file in a submitted archive must be accompanied by its SHA1 **uncompressed** content digest file, digitally signed with the submitter's signing key.

```
# Signed digest of uncompressed contents, will be dlat-p--001-20120920--mycode.txt.shal
openssl dgst -shal -sign subKey.pem -out dlat-p--001-20120920--mycode.txt.shal \
dlat-p--001-20120920--mycode.txt

# Signed digest of already compressed contents
bunzip2 -c dlat-p--001-20120920--mycode.txt.bz2 | openssl dgst -shal -sign subKey.pem \
-out dlat-p--001-20120920--mycode.txt.shal

# Verification of uncompressed contents using
# signed digest dlat-p--001-20120920--mycode.txt.shal
openssl dgst -shal -verify subKey.pem.pub -signature dlat-p--001-20120920--mycode.txt.shal \
dlat-p--001-20120920--mycode.txt

# Verification of compressed contents
bunzip2 -c dlat-p--001-20120920--mycode.txt.bz2 | openssl dgst -shal -verify subKey.pem.pub \
-signature dlat-p--001-20120920--mycode.txt.shal
```

The procedure to submit analyzed contents to BLUEPRINT DCC also involves first having the raw data used for the analysis in the European Genome-phenome Archive (EGA), as all the metadata entries from the analyzed contents to be stored in BLUEPRINT DCC must point to the original raw data.

1.2.1 File Naming Conventions

Submitted files, containing either sample tracking data or analyzed experiment contents, must follow next file naming convention

```
featureType-fileType--institutionCode-dateFileCreated--freeField.txt
```

```
featureType-fileType--institutionCode-dateFileCreated--freeField.txt.sha1
```

The file name components are mapped in the next way:

Components	Description	Key
	Sample Tracking data	sdata
	Gene Expression	exp
fo atura Tuna	Exon Junctions	jcn
featureType	DNA *-lation (Methylation, Hydroxy-Methylation,	dlat
	Formylation, etc)	
	Protein-DNA interactions	pdna
	Regulatory regions	rreg
	Metadata file	m
	Primary data file	р
	Secondary data file	S
filoTupo	Gene expression file	g
fileType	Donor file	donor
	Specimen file	specimen
	Sample file	sample
	Donor's Family file	family
institutionCode	Institution submitting data	CV Table A.5.1
dateFileCreated	The date on which the file is created	YYYYMMDD (ISO-8601)
freeField	An alphanumeric field (max length of 16 characters)	e.g.: mysample, 0B1845J
	where submitters can put internal codes, file sequence	
	numbers, etc	

Different file types of the same feature type are interrelated, because the data they are storing is intertwined. Specific relations are defined on the documentation of each feature type and their file types. For instance, information stored in a primary data file is related and depends on the data from its corresponding metadata file, and the same happens to secondary data files and primary data files. Metadata file contents are related to sample tracking data sample files.

1.2.2 Tabular File Structure

The submitted analyzed contents are kept in tab-delimited text files. General comments may be added to the beginning of the file with a hash ('#') prefixed at beginning of each comment line. The first non-comment line is the header containing the names of the columns. Each column corresponds to a data element defined in DCC Submission Tabular Formats specification (Chapter 2).

There is a subset of comment lines used to attach data labels to the text files. These data labels follow the form '##labelName value [value ...]'. Currently acknowledged data labels are:

- **format**: This label is **required**, and its value defines the BLUEPRINT data formatting schema used on the file.
- depends: Although this label is not always required, it is important to validate the data coherence of the whole data set, because it ensures related data is not corrupted. The values of this label are the file on the same submission this file is related to (for instance, the name of a metadata file), and the SHA1 digest value (in its hexadecimal representation) of that file's contents.

There are several ways to generate the SHA1 digest of a file, like libraries in most of the programming languages and command-line tools:

```
# Getting the SHA1 digest value of uncompressed contents using OpenSSL
openssl dgst -sha1 dlat-p--001-20120920--mycode.txt

SHA1(dlat-p--001-20120920--mycode.txt)= 81ae49a7014d2d0260625d3535fa6e2a4a0bc06f

# Getting the SHA1 digest value of uncompressed contents using sha1sum
sha1sum dlat-p--001-20120920--mycode.txt

81ae49a7014d2d0260625d3535fa6e2a4a0bc06f dlat-p--001-20120920--mycode.txt
```

An example file is shown below (note that parts of the lines are omitted for readability):

```
# This is an example of a primary analysis file for simple somatic mutations.
# File name: dlat-p--001-20120920--mycode.txt
#
# And it has its labels
##format 0.2.3.2-566596e769717c140078683abfe46f007f199452
##depends dlat-m--001-20120920--mycode.txt 03366af5145107cc818f4827e86b61dcf998ff29
            ⇒analyzed_sample_id
                                ⇒d_lated_fragment_id
                                                    ⇒chromosome
                                                                 \rightarrow...
                                                                        ⊣note
analvsis_id
→dlat:001:1234ff33
                                                    ∌1
                                                           \rightarrow...
                                                                  ⇒#FF#
→dlat:001:00019878
                                                    ∃1
                                                           \dashv...
                                                                  →#FF#
                                                    ⇒|21
                                                           ⇒...
⇒dlat:001:a712838
                                                                  →#FF#
⇒dlat:001:abebdZZZZZ
                                                                  ⇒#FF#
```

All the declared columns for each file type must be set. Data columns are labeled as identifier or reference (I), required (R), desirable (D) or optional (O). Data providers (i.e. submitters) must put all the efforts in order to provide values for the idref and required data columns. The exception for this rule are the desirable fields, required fields which can be unknown on the first submissions, but in that case the fields these exceptions are properly documented.

There are several possible reasons why a column value (either desirable or optional) has not been provided. Next reserved codes must be used to describe the reason:

Code	Meaning
#FF#	Data not supplied at this time (for future fill)
#NA#	Not applicable for the context of the surrounding knowledge
#V0#	Data verified to be unknown (void, undef, null)
#DE#	Data derived from a required or idref field

Some data columns described in this submission manual contain values used as identifiers on BLUEPRINT DCC (e.g. analysis_id, regulatory_region_id, ...). As such, these identifiers should uniquely identify the entity they are referring (an analysis, a regulatory region, ...), and the identifier's value should be globally unique within a center's data submission. Also, these identifiers should be consistent along the different data submissions and releases. If you have to generate your own identifiers, there are some general recommendations, like using the same prefix for the identifiers of the same kind.

When you are submitting string values for columns which can contain URLs or multiple values delimited by commas, each separate value string, before being joined, should be URI encoded.

1.3 File Submission Procedure

Files with the contents to be submitted, along with their corresponding signed disgest, must be sent in a single tar archive. Either the tar archive or its embedded contents should be submitted compressed, using gzip, bzip2 or xz formats.

To be finished/defined

Chapter 2

DCC Submission Tabular Formats

2.1 Sample Tracking Submission File Specifications

Overview

There are three required sample and tissue annotation submission files, and one optional template file:

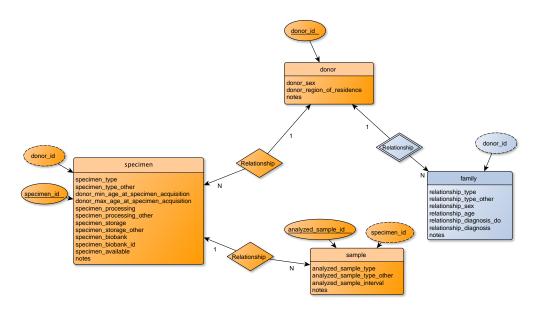


Figure 2.1: Sample Tracking Data Sub-Schema

Core Sample Tracking Data Files

- 1. Donor Data File (donor)

 Mandatory information about the donor's age, gender and vital status.
- 2. Specimen Data File (specimen)

 Mandatory information about a specimen that was obtained from a donor. There may be several specimens per donor that were obtained concurrently or at different times.
- 3. Analyzed Sample Data File (sample)

 Mandatory information about an analyzed sample that was subjected to molecular analysis. There may be several analyzed samples per specimen, for example, blood samples at different ages.

All data submissions to the DCC must include the three core sample tracking data files.

Optional Template Files

Donor Family History (family)
 Optional details about family history of the donor

Coding of donor IDs

The three mandatory data files contain donor, specimen and analyzed sample IDs, respectively. These IDs are to be coded specifically for BLUEPRINT purposes and only the submitting group will keep the key that will permit to link back the data to the individual donors. The key must not be communicated to the data users. It should not be derived from other IDs such as biobank or hospital identifiers. These IDs are to be coded in such a way that they cannot be tracked back to the individual donors, except by the submitting group. IDs are assigned by each submitting group, and must be unique within all the data submitted by that group (i.e. no duplicate IDs allowed). The DCC will prevent collisions between similar IDs submitted by different groups by including the project source column by default in all BioMart queries.

2.1.1 Donor Data File

Donor Data File [donor] (required)

This submission file describes a donor from which one or more specimens were obtained.

Table 2.1.1: Donor Data File

Name	Туре	Need	Description / Values
donor_id	string	I	Unique identifier for the donor; assigned by data provider.
donor_region_of_residence	<pre>string[] (array seps ,)</pre>	R	Country, and optionally state or province code, but not city. ISO3166-1-alpha-2 or ISO3166-2 codes, eg: "CA" or "CA-ON" (See external CV description A.4)
donor_sex	integer	D	Donor biological sex. "Other" has been removed from the controlled vocabulary due to identifiability concerns. 1 = male 2 = female
notes	string	0	Any additional non-identifying information can be included here.

2.1.2 Specimen Data File

Specimen Data File [specimen] (required)

This submission file describes a specimen from which one or more samples were derived. Use additional rows for more than one specimen from the same patient. If more than one specimen was extracted during the same procedure, each gets a distinct ID.

Table 2.1.2: Specimen Data File

Name	Type	Need	Description / Values
specimen_id	string	I	Unique identifier for the specimen assigned by data provider.
			Continued on next page

Table 2.1.2 – continued from previous page

Table 2.1.2 – continued from previous page							
Name	Туре	Need	Description / Values				
donor_id	string	Ι	Unique identifier for the donor; assigned by data provider.				
donor_max_age_at_specimen_acquisition	duration	D	Donor max age when the specimen was acquired. If it is '#DE#', then it is the same value as "donor_min_age_at_specimen_acquisition"				
donor_min_age_at_specimen_acquisition	duration	R	Donor minimal age when the specimen was acquired, in ISO-8601 duration (basic format)				
notes	string	0	Any additional non-identifying information can be included here.				
specimen_available	boolean	0	Whether additional tissue is available for followup studies.				
specimen_biobank	string	0	If the specimen was obtained from a biobank, provide the biobank name here				
specimen_biobank_id	string	0	If the specimen was obtained from a biobank, provide the biobank accession number here.				
specimen_processing	string	R	Description of technique used to process specimen 1 = cryopreservation in liquid nitrogen (dead tissue) 2 = cryopreservation in dry ice (dead tissue) 3 = cryopreservation of live cells in liquid nitrogen 4 = cryopreservation, other 5 = formalin fixed, unbuffered 6 = formalin fixed, buffered 7 = formalin fixed & paraffin embedded 8 = fresh 9 = other technique				
specimen_processing_other	string	0	If "other" specified for specimen_processing, may indicate technique here.				
specimen_storage	string	R	Description of how specimen was stored. For specimens that were extracted freshly or immediately cultured, answer (1) "NA". 1 = frozen, liquid nitrogen 2 = frozen, -70 freezer 3 = frozen, vapor phase 4 = RNA later frozen 5 = paraffin block 6 = cut slide 7 = other				
			Continued on next page				

Table 2.1.2 – concluded from previous page

Name	Туре	Need	Description / Values
specimen_storage_other	string	0	If "other" specified for specimen_storage, may indicate technique here.
specimen_type	integer	R	Controlled vocabulary description of specimen type. 1 = primary tumour 2 = tumour local recurrence 3 = tumour metastasis to local lymph node 4 = tumour metastasis to distant location 5 = peripheral blood 6 = bone marrow 7 = lymph node 8 = normal control (tissue adjacent to primary) 9 = normal control (blood) 10 = normal control (other) 11 = disease tissue (other) 12 = cord blood
specimen_type_other	string	0	Free text description of site of specimen if "normal control (other)" or "disease tissue (other)" was specified in specimen_type field.

2.1.3 Analyzed Sample Data File

Analyzed Sample Data File [sample] (required)

This submission file describes an analyzed sample on which molecular characterization was performed. It includes both control samples (from healthy people) and samples from ill people.

Table 2.1.3: Analyzed Sample Data File

Name	Туре	Need	Description / Values
analyzed_sample_id	string	I	Unique identifier for the sample assigned by data provider
specimen_id	string	I	Unique identifier for the specimen assigned by data provider.
analyzed_sample_interval	integer	0	Interval from specimen acquisition to sample use in an analytic procedure (e.g. DNA extraction), in days
analyzed_sample_type	integer	R	Controlled vocabulary description of sample type 1 = Normal blood 2 = Leukemic blood 3 = Normal control adjacent to primary 4 = Normal control from non-tumour site 5 = Control from cell line derived from normal tissue 6 = Normal mouse host 7 = Primary tumour 8 = Mouse xenograft derived from tumour 9 = Cell line derived from tumour 10 = Cell line derived from xenograft 11 = Other (specify)
analyzed_sample_type_other	string	0	Free text description of site of sample if "other" was specified in sample_type field
notes	string	0	Any additional non-identifying information can be included here.
purified_cell_type	string	R	Purified cell type for the sample (See external CV description A.3)

2.2 Gene Expression

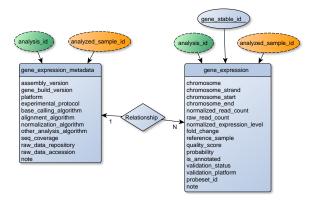


Figure 2.2: Gene Expression Sub-Schema

2.2.1 Expression - Metadata File

Expression [exp] – Metadata File [m]

Table 2.2.1: Expression - Metadata File

Table 2.2.1: Expression - Metadata File							
Name	Туре	Need	Description / Values				
analysis_id	string	Ι	Unique identifier for the analysis performed for a particular set of samples				
analyzed_sample_id	string	I	Unique identifier for the analyzed sample				
alignment_algorithm	compound name;url	R	Name of alignment algorithm and URL to written protocol				
analysis_group_id	integer	R	Identifier of the analysis group (i.e. the one who pre- pared/run the pipeline) (See CV A.5)				
assembly_version	integer	R	Version of reference genome assembly (See CV A.8)				
data_status	integer	R	The status of the analysis over the associated raw data 0 = Raw data available at the EGA, but not more 1 = Raw data available at the EGA, analysis in process 2 = Analysis results obtained (analysis finished)				
experimental_group_id	integer	R	Identifier of the experimental group who did the experimental analysis (See CV A.5)				
experimental_protocol	compound name;url	0	Name of experimental protocol and URL to written protocol				
note	string	0	Optional field to leave notes				
other_analysis_algorithm	<pre>compound[] name;url (array seps ,)</pre>	0	Names of other analysis algorithms. Separate multiple algorithms by commas.				
platform	integer	R	Platform or technology used in the detection phase (See CV A.6)				
raw_data_accession	<pre>compound accession;url</pre>	0	Accession and URL for referencing the raw data at the public repository				
raw_data_repository	integer	R	Public repository where raw data is submitted (#) (See CV A.9)				
seq_coverage	decimal	0	Sequence coverage if analyzed by sequencing plat- forms				

2.2.2 Expression - Gene File

 $Expression \ [exp]-Gene \ File \ [g]$

Table 2.2.2: Expression - Gene File

Name	Туре	Need	Description / Val	ues
		(Continued on next p	age

Table 2.2.2 – continued from previous page


	Table 2.2.2 – continued from previous page					
Name	Туре	Need	Description / Values			
analysis_id	string	I	Unique identifier for the analysis performed for a particular set of samples			
analyzed_sample_id	string	I	Unique identifier for the analyzed sample			
gene_stable_id	string	I	For annotated gene, use Ensembl gene ID. Otherwise, use assemblyBuild_chr_start_end where assemblyBuild is hg18 or hg19.			
chromosome	integer	R	Name of the chromosome containing the experimentally detected feature (mutation, variation, expression,) (See CV A.7)			
chromosome_end	integer	R	End position of the mutation/variation on the chromosome			
chromosome_start	integer	R	Start position of the mutation/variation on the chromosome			
chromosome_strand	integer	0	Strand where it was detected the mutation/variation on the chromosome -1 = Reverse strand 1 = Forward strand			
fold_change	decimal	0	Expressed fold change if differential expression is measured			
is_annotated	boolean	R	If it is true, it indicate if the expressed fragment is annotated in Ensembl (i.e. gene_stable_id contains a Ensembl Gene Identifier)			
normalized_expression_level	decimal	0	Normalized value of expression level if analyzed by microarray platforms			
normalized_read_count	decimal	R	Normalized count of sequencing reads if analyzed by sequencing platforms			
note	string	0	Optional field to leave notes			
probability	decimal	0	Probability of the mutation/variation call			
probeset_id	string	0	ID of the probeset used in microarray if analyzed by microarray platform			
quality_score	decimal	0	Average quality score for the mutation/variation call			
raw_read_count	integer	R	Raw count of sequencing reads if analyzed by sequencing platforms			
	1	1	Continued on next page			

Table 2.2.2 – concluded from previous page

Name	Туре	Need	Description / Values
reference_sample	string	0	ID of the reference analyzed sample if differential expression is measured
validation_platform	integer	0	Platform or technology used in validation (See CV A.6)
validation_status	integer	R	Indicate if the mutation/variation has been validated -1 = Not valid 0 = Not tested 1 = Validated

2.3 Exon Junction

The following diagram, based on the one from ICGC DCC manual, illustrates how junction_id should be generated, how junction_read_count, exon1_number_bases and exon2_number_bases are calculated:

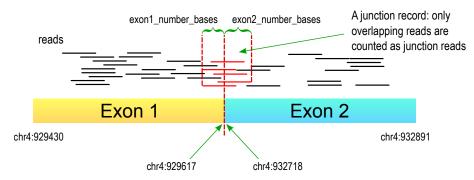

junction_id would be: hg19_4_929617_4_932718 junction read count is: 5

Figure 2.3: Junction Read Count explanation

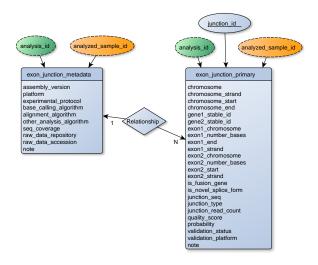


Figure 2.4: Exon Junction Sub-Schema

2.3.1 Exon Junction - Metadata File

Exon Junction [jcn] – Metadata File [m]

Table 2.3.1: Exon Junction - Metadata File

Name	Tupe	Need	ion - Metadata File Description / Values
analysis_id	string	I	Unique identifier for the analysis performed for a par-
	Jei Ing	-	ticular set of samples
analyzed_sample_id	string	I	Unique identifier for the analyzed sample
alignment_algorithm	compound name;url	R	Name of alignment algorithm and URL to written protocol
analysis_group_id	integer	R	Identifier of the analysis group (i.e. the one who pre- pared/run the pipeline) (See CV A.5)
assembly_version	integer	R	Version of reference genome assembly (See CV A.8)
data_status	integer	R	 The status of the analysis over the associated raw data 0 = Raw data available at the EGA, but not more 1 = Raw data available at the EGA, analysis in process 2 = Analysis results obtained (analysis finished)
experimental_group_id	integer	R	Identifier of the experimental group who did the experimental analysis (See CV A.5)
experimental_protocol	compound name;url	0	Name of experimental protocol and URL to written protocol
note	string	0	Optional field to leave notes
other_analysis_algorithm	<pre>compound[] name;url (array seps ,)</pre>	0	Names of other analysis algorithms. Separate multiple algorithms by commas.
platform	integer	R	Platform or technology used in the detection phase (See CV A.6)
raw_data_accession	<pre>compound accession;url</pre>	0	Accession and URL for referencing the raw data at the public repository
raw_data_repository	integer	R	Public repository where raw data is submitted (#) (See CV A.9)
seq_coverage	decimal	0	Sequence coverage if analyzed by sequencing plat- forms
		<u>I</u>	Continued on next page

2.3.2 Exon Junction - Primary Analysis File

Exon Junction [jcn] – Primary Analysis File [p]

Table 2.3.2: Exon Junction - Primary Analysis File

Table 2.3.2: Exon Junction - Primary Analysis File					
Name	Type	Need	Description / Values		
analysis_id	string	I	Unique identifier for the analysis performed for a particular set of samples		
analyzed_sample_id	string	I	Unique identifier for the analyzed sample		
junction_id	string	I	For known exons, use exonID1_exonID2 where exonID1 and exonID2 are Ensembl IDs of the 5' and 3' exons, respectively. For novel or putative exons, use assembly—Build_exon1chr_exon1end_exon2chr_exon2start where assemblyBuild is hg18 or hg19; exon1chr and exon2chr are the chromosomes of the 5' and 3' exons, respectively; exon1end is the end position of the 5' exon; exon2start is the start position of the 3' exon.		
chromosome	integer	R	Name of the chromosome containing the experimentally detected feature (mutation, variation, expression,) (See CV A.7)		
chromosome_end	integer	R	End position of the mutation/variation on the chromosome		
chromosome_start	integer	R	Start position of the mutation/variation on the chromosome		
chromosome_strand	integer	0	Strand where it was detected the mutation/variation on the chromosome -1 = Reverse strand 1 = Forward strand		
exon1_chromosome	integer	R	Name of the chromosome containing the 5' exon (#) (See CV A.7)		
exon1_end	integer	R	End position of the 5' exon on the chromosome		
exon1_number_bases	integer	R	Number of bases from 5' exon		
exon1_strand	integer	0	Chromosome strand of the 5' exon -1 = Reverse strand 1 = Forward strand Continued on next page		

Table 2.3.2 – continued from previous page

Name	Type	2.3.2 – 0 Need	ontinued from previous page Description / Values
	· ·		
exon2_chromosome	integer	R	Name of the chromosome containing the 3' exon (#) (See CV A.7)
exon2_number_bases	integer	R	Number of bases from 3' exon
e.conz_number_buses	Integer	١٠,	Trainer of buses from 5 exon
exon2_start	integer	R	Start position of the 3' exon on the chromosome
exon2_strand	integer	0	Chromsome strand of the 3' exon
			-1 = Reverse strand
			1 = Forward strand
gene1_stable_id	string	R	Stable ID of the gene containing the 5' exon at
			the junction. For annotated gene, use Ensembl
			gene ID. For putative and novel gene, use assem-
			blyBuild_chr_start_end where assemblyBuild can be
			hg18 or hg19.
gana? atal-!-!		0	In the case of a fusion game arrest to the Ctable ID of
gene2_stable_id	string	0	In the case of a fusion gene, provide the Stable ID of
			the gene containing the 3' exon at the junction. For annotated genes, use Ensembl gene ID. For putative
			annotated genes, use Ensembl gene ID. For putative and novel genes, use assemblyBuild_chr_start_end
			where assemblyBuild can be hg18 or hg19.
			where assembly buttu can be higher or higher.
is_fusion_gene	boolean	0	Indicate if the function is the result of a fusion gene
is_rusion_gene	Doocean	J	mateure if the function is the result of a fusion gene
is_novel_splice_form	boolean	0	Indicate if the splice form is novel
→ - · -			'
junction_read_count	integer	R	Count of sequencing reads that span across exons
junction_seq	string	0	Provide junction sequence if either is_fusion_gene or
			is_novel_splice_form is true
junction_type	integer	0	Type of junction
			1 = Canonical
			2 = Non-canonical
		-	3 = U12
note	string	0	Optional field to leave notes
	42-7		Duck abilian of the manager of the coll
probability	decimal	0	Probability of the mutation/variation call
quality_score	decimal	0	Average quality score for the mutation/variation call
quality_score	uecimat	U	Average quality score for the mutation/variation call
validation_platform	integer	0	Platform or technology used in validation
vattaatton_ptationiii	In cage	J	(See CV A.6)
validation_status	integer	R	Indicate if the mutation/variation has been validated
validation_status	In cage	١١.	-1 = Not valid
			0 = Not tested
			1 = Validated
	l .	<u> </u>	Continued on next page
			. 1 3 .

Name Type Need Description / Values

2.4 Protein-DNA interactions

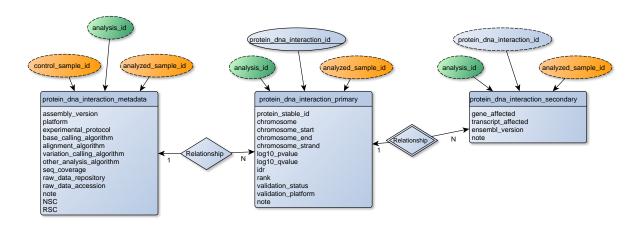


Figure 2.5: Protein-DNA interactions Sub-Schema

2.4.1 Protein-DNA interaction - Metadata File

Protein-DNA [pdna] - Metadata File [m]

Table 2.4.1: Protein-DNA interaction - Metadata File					
Name	Туре	Need	Description / Values		
analysis_id	string	I	Unique identifier for the analysis performed for a particular set of samples		
analyzed_sample_id	string	I	Unique identifier for the analyzed sample		
NSC	decimal	0	Normalized strand-cross correlation of the analysis		
RSC	decimal	0	Relative strand-cross correlation of the analysis		
alignment_algorithm	compound name;url	R	Name of alignment algorithm and URL to written protocol		
analysis_group_id	integer	R	Identifier of the analysis group (i.e. the one who pre- pared/run the pipeline) (See CV A.5)		
assembly_version	integer	R	Version of reference genome assembly (See CV A.8)		
			Continued on next page		

Table 2.4.1 – concluded from previous page

Name	Туре	Need	Description / Values
control_sample_id	string	R	Unique identifier for the analyzed control/matched sample
data_status	integer	R	The status of the analysis over the associated raw data 0 = Raw data available at the EGA, but not more 1 = Raw data available at the EGA, analysis in process 2 = Analysis results obtained (analysis finished)
experimental_group_id	integer	R	Identifier of the experimental group who did the experimental analysis (See CV A.5)
experimental_protocol	compound name;url	0	Name of experimental protocol and URL to written protocol
note	string	0	Optional field to leave notes
other_analysis_algorithm	<pre>compound[] name;url (array seps ,)</pre>	0	Names of other analysis algorithms. Separate multiple algorithms by commas.
platform	integer	R	Platform or technology used in the detection phase (See CV A.6)
raw_data_accession	<pre>compound accession;url</pre>	0	Accession and URL for referencing the raw data at the public repository
raw_data_repository	integer	R	Public repository where raw data is submitted (#) (See CV A.9)
seq_coverage	decimal	0	Sequence coverage if analyzed by sequencing plat- forms

2.4.2 Protein-DNA interaction - Primary Analysis File

 $Protein-DNA \left[pdna\right]-Primary \ Analysis \ File \left[p\right]$

Table 2.4.2: Protein-DNA interaction - Primary Analysis File

Name	Туре	Need	Description / Values
analysis_id	string	I	Unique identifier for the analysis performed for a particular set of samples
analyzed_sample_id	string	I	Unique identifier for the analyzed sample
protein_dna_interaction_id	string	I	Unique identifier for the protein-DNA interaction
Continued on next pag			

Table 2.4.2 – concluded from previous page

Name	Type	Need	Description / Values
chromosome	integer	R	Name of the chromosome containing the experimentally detected feature (mutation, variation, expression,) (See CV A.7)
chromosome_end	integer	R	End position of the mutation/variation on the chromosome
chromosome_start	integer	R	Start position of the mutation/variation on the chromosome
chromosome_strand	integer	0	Strand where it was detected the mutation/variation on the chromosome -1 = Reverse strand 1 = Forward strand
idr	decimal	D	Irreproducible discovery rate
log10_pvalue	decimal	R	-log10(p-value)
log10_qvalue	decimal	D	-log10(q-value) , which available for peaks, but not for broad peaks
note	string	0	Optional field to leave notes
protein_stable_id	string	R	Stable id of the interacting protein, antibody or protein complex
rank	<pre>compound[] rank;value (array seps)</pre>	0	Kind of used ranking and its value, in the form "rank;value". As it can hold more than one value, they are separated by bars
validation_platform	integer	0	Platform or technology used in validation (See CV A.6)
validation_status	integer	R	Indicate if the mutation/variation has been validated -1 = Not valid 0 = Not tested 1 = Validated

2.4.3 Protein-DNA interaction - Secondary Analysis File

 $Protein-DNA \ [pdna] - Secondary \ Analysis \ File \ [s]$

Table 2.4.3: Protein-DNA interaction - Secondary Analysis File

Name	Type	Need	Description / Values
analysis_id	string	I	Unique identifier for the analysis performed for a particular set of samples
			Continued on next page

Table 2.4.3 – concluded from previous page

Name	Type	Need	Description / Values
analyzed_sample_id	string	I	Unique identifier for the analyzed sample
ensembl_version	integer	R	Version of Ensembl gene build used for annotation
gene_affected	<pre>string[] (array seps)</pre>	R	Gene affected. Use Ensembl gene id, separated by when there is more than one. If no gene is affected, don't put an entry (See external CV description A.1)
note	string	0	Optional field to leave notes
protein_dna_interaction_id	string	I	Unique identifier for the protein-DNA interaction
transcript_affected	<pre>string[] (array seps)</pre>	0	Transcript on the protein-DNA interaction area. Use Ensembl transcript id. Separate multiple transcripts with vertical bars in the form of transcriptA transcriptB transcriptC (See external CV description A.2)

2.5 Regulatory Regions

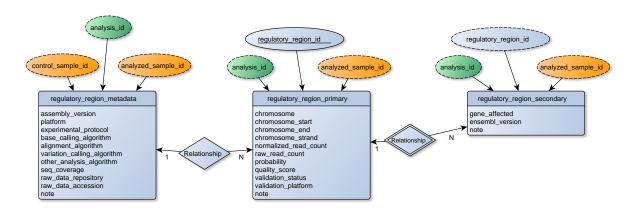


Figure 2.6: Regulatory Regions Sub-Schema

2.5.1 Regulatory regions - Metadata File

Regulatory regions [rreg] – Metadata File [m]

Table 2.5.1: Regulatory regions - Metadata File

Name	Type	Need	Description / Values
analysis_id	string	I	Unique identifier for the analysis performed for a particular set of samples
			Continued on next page

Table 2.5.1 – concluded from previous page

Name	Type	Need	Description / Values
analyzed_sample_id	string	I	Unique identifier for the analyzed sample
unutgzeu_sumpte_tu	Stiring	-	onique tuentater for the unatyzed sumpte
alignment_algorithm	compound	R	Name of alignment algorithm and URL to written pro-
	name;url		tocol
analysis_group_id	integer	R	Identifier of the analysis group (i.e. the one who pre-
			pared/run the pipeline)
			(See CV A.5)
assembly_version	integer	R	Version of reference genome assembly
	_		(See CV A.8)
control_sample_id	string	R	Unique identifier for the analyzed control/matched
			sample
data_status	intonon	R	The status of the analysis over the associated raw data
udta_status	integer	K	0 = Raw data available at the EGA, but not
			more
			1 = Raw data available at the EGA, analysis in
			process
			2 = Analysis results obtained (analysis fin-
			ished)
experimental_group_id	integer	R	Identifier of the experimental group who did the ex-
			perimental analysis
			(See CV A.5)
experimental_protocol	compound	0	Name of experimental protocol and URL to written
	name;url		protocol
noto		0	Ontional field to leave notes
note	string	0	Optional field to leave notes
other_analysis_algorithm	compound[]	0	Names of other analysis algorithms. Separate multi-
zaatgoto_atgottiiiii	name;url		ple algorithms by commas.
	(array seps)		. ,
platform	integer	R	Platform or technology used in the detection phase
			(See CV A.6)
raw_data_accession	compound	0	Accession and URL for referencing the raw data at the
	accession;url		public repository
			B.H.
raw_data_repository	integer	R	Public repository where raw data is submitted (#)
000 0000000	da a 2 1	0	(See CV A.9)
seq_coverage	decimal	0	Sequence coverage if analyzed by sequencing plat- forms
			IUTIIIS

2.5.2 Regulatory regions - Primary Analysis File

 $Regulatory\ regions\ [rreg]-Primary\ Analysis\ File\ [p]$

Table 2.5.2: Regulatory regions - Primary Analysis File

Name Type Need Description / Values

Continued on next page

Table 2.5.2 – concluded from previous page

Name Type Need Description / Values				
	Туре	- Neeu	<u> </u>	
analysis_id	string	I	Unique identifier for the analysis performed for a particular set of samples	
analyzed_sample_id	string	I	Unique identifier for the analyzed sample	
regulatory_region_id	string	I	Unique identifier for the identified regulatory region	
chromosome	integer	R	Name of the chromosome containing the experimentally detected feature (mutation, variation, expression,) (See CV A.7)	
chromosome_end	integer	R	End position of the mutation/variation on the chromosome	
chromosome_start	integer	R	Start position of the mutation/variation on the chromosome	
chromosome_strand	integer	0	Strand where it was detected the mutation/variation on the chromosome -1 = Reverse strand 1 = Forward strand	
normalized_read_count	decimal	R	Normalized count of sequencing reads if analyzed by sequencing platforms	
note	string	0	Optional field to leave notes	
probability	decimal	0	Probability of the mutation/variation call	
quality_score	decimal	0	Average quality score for the mutation/variation call	
raw_read_count	integer	R	Raw count of sequencing reads if analyzed by sequencing platforms	
validation_platform	integer	0	Platform or technology used in validation (See CV A.6)	
validation_status	integer	R	Indicate if the mutation/variation has been validated -1 = Not valid 0 = Not tested 1 = Validated	

2.5.3 Regulatory regions - Secondary Analysis File

 $Regulatory\ regions\ [rreg]-Secondary\ Analysis\ File\ [s]$

Table 2.5.3: Regulatory regions - Secondary Analysis File

Name	Туре	Need	Description / Values
analysis_id	string	I	Unique identifier for the analysis performed for a particular set of samples
			Continued on next page

Table 2.5.3 – concluded from previous page

Name	Туре	Need	Description / Values
analyzed_sample_id	string	I	Unique identifier for the analyzed sample
ensembl_version	integer	R	Version of Ensembl gene build used for annotation
gene_affected	<pre>string[] (array seps)</pre>	R	Gene affected. Use Ensembl gene id, separated by when there is more than one. If no gene is affected, don't put an entry (See external CV description A.1)
note	string	0	Optional field to leave notes
regulatory_region_id	string	I	Unique identifier for the identified regulatory region

2.6 DNA *-lation (Methylation, Hydroxy-Methylation, Formylation, etc...)

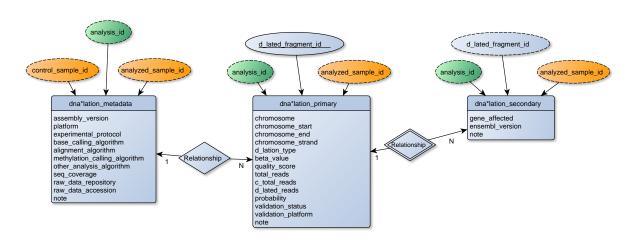


Figure 2.7: DNA Methylation , Hydroxy-Methylation, Formylation, etc... Sub-Schema

2.6.1 DNA *-lation - Metadata File

DNA *-lation [dlat] - Metadata File [m]

Name	Туре	Need	Description / Values
analysis_id	string	I	Unique identifier for the analysis performed for a particular set of samples
analyzed_sample_id	string	Ι	Unique identifier for the analyzed sample
Continued on next page			

Table 2.6.1 – concluded from previous page

Name Type Need Description / Values						
	Туре		·			
alignment_algorithm	compound	R	Name of alignment algorithm and URL to written pro-			
	name;url		tocol			
		_				
analysis_group_id	integer	R	Identifier of the analysis group (i.e. the one who pre-			
			pared/run the pipeline)			
			(See CV A.5)			
assembly_version	integer	R	Version of reference genome assembly			
		-	(See CV A.8)			
control_sample_id	string	R	Unique identifier for the analyzed control/matched sample			
			Sample			
data_status	integer	R	The status of the analysis over the associated raw data			
data_statas	integer	.,	0 =Raw data available at the EGA, but not			
			more 1 = Raw data available at the EGA, analysis in			
			process			
			2 = Analysis results obtained (analysis fin-			
			ished)			
experimental_group_id	integer	R	Identifier of the experimental group who did the ex-			
			perimental analysis			
			(See CV A.5)			
experimental_protocol	compound	0	Name of experimental protocol and URL to written			
	name;url		protocol			
methylation_calling_algorithm	compound	0	Name of variation calling algorithm and URL to writ-			
	name;url		ten protocol			
		_				
note	string	0	Optional field to leave notes			
othor analysis alassithus		0	Names of other analysis also with as Consust will			
other_analysis_algorithm	compound[]	0	Names of other analysis algorithms. Separate multi-			
	<pre>name;url (array seps ,)</pre>		ple algorithms by commas.			
platform	integer	R	Platform or technology used in the detection phase			
padom	Integer	IX.	(See CV A.6)			
raw_data_accession	compound	0	Accession and URL for referencing the raw data at the			
	accession;url	_	public repository			
	· • · ·		, , ,			
raw_data_repository	integer	R	Public repository where raw data is submitted (#)			
. 3	-		(See CV A.9)			
seq_coverage	decimal	0	Sequence coverage if analyzed by sequencing plat-			
			forms			

2.6.2 DNA *-lation - Primary Analysis File

DNA *-lation [dlat] - Primary Analysis File [p]

Table 2.6.2: DNA *-lation - Primary Analysis File

Name Type Need Description / Values

Continued on next page

Table 2.6.2 – continued from previous page

Table 2.6.2 – continued from previous page					
Name	Туре	Need	Description / Values		
analysis_id	string	I	Unique identifier for the analysis performed for a particular set of samples		
analyzed_sample_id	string	I	Unique identifier for the analyzed sample		
d_lated_fragment_id	string	I	Unique identifier for the methy- lated fragment, in the form d'*lationType chromosome_chromosomeStart_chromosome		
beta_value	decimal	0	DNA *-lation beta value for interrogated site		
c_total_reads	decimal	R	Reads which has identified this position as a cytosine		
chromosome	integer	R	Name of the chromosome containing the experimentally detected feature (mutation, variation, expression,) (See CV A.7)		
chromosome_end	integer	R	End position of the mutation/variation on the chromosome		
chromosome_start	integer	R	Start position of the mutation/variation on the chromosome		
chromosome_strand	integer	0	Strand where it was detected the mutation/variation on the chromosome -1 = Reverse strand 1 = Forward strand		
d_lated_reads	decimal	R	Reads which has identified this position as a DNA *lated cytosine		
d_lation_type	string	R	Type of DNA *-lation (Methylation, Hydroxy-Methylation, Formylation, etc) m = Methylation (cytosine) hm = Hydroxy-Methylation (cytosine) hmU = Hydroxy-Methylation (uracil) f = Formylation (cytosine) ca = Carboxylation (cytosine)		
note	string	0	Optional field to leave notes		
probability	decimal	0	Probability of the DNA *-lation call		
quality_score	decimal	0	Quality score for the DNA *-lation call		
total_reads	decimal	R	Total number of reads over this position/segment, either identifying or not a cytosine, for sequencing platforms. Mean reads depth for other technologies		

Table 2.6.2 – concluded from previous page

Name	Туре	Need	Description / Values
validation_platform	integer	0	Platform or technology used in validation (See CV A.6)
validation_status	integer	R	Indicate if the mutation/variation has been validated -1 = Not valid 0 = Not tested 1 = Validated

2.6.3 DNA *-lation - Secondary Analysis File

DNA *-lation [dlat] - Secondary Analysis File [s]

Table 2.6.3: DNA *-lation - Secondary Analysis File

Name	Type	Need	Description / Values	
analysis_id	string	I	Unique identifier for the analysis performed for a particular set of samples	
analyzed_sample_id	string	I	Unique identifier for the analyzed sample	
d_lated_fragment_id	string	I	Unique identifier for the methy- lated fragment, in the form d'*lationType chromosome_chromosomeStart_chromosom	ieEnd
ensembl_version	integer	R	Version of Ensembl gene build used for annotation	
gene_affected	<pre>string[] (array seps)</pre>	R	Gene affected. Use Ensembl gene id, separated by when there is more than one. If no gene is affected, don't put an entry (See external CV description A.1)	
note	string	0	Optional field to leave notes	

Appendix A

Controlled Vocabularies

A.1 Ensembl Genes

Valid Ensembl Genes identifiers
(See it at http://jan2013.archive.ensembl.org/Homo_sapiens/Info/Index)

A.2 Ensembl Transcripts

Valid Ensembl Transcript identifiers
(See it at http://jan2013.archive.ensembl.org/Homo_sapiens/Info/Index)

A.3 Cell Ontology

The Cell Ontology is designed as a structured controlled vocabulary for cell types (See it at http://cellontology.org/)

A.4 ISO 3166-1 and ISO 3166-2

ISO 3166 is the International Standard for country codes and codes for their subdivisions. The purpose of ISO 3166 is to establish internationally recognised codes for the representation of names of countries, territories or areas of geographical interest, and their subdivisions.

(See them at http://www.iso.org/iso/country_codes.htm and http://en.wikipedia.org/wiki/ISO_ 3166-2)

Alias	Key	: ISO 3166-1 and ISO 3166-2 aliases Description
ALIAS:EAL	GB-CAM GB-ESS GB-HRT GB-NFK GB-SFK	East Anglia: United Kingdom region composed of the administrative counties of Norfolk to the north, Suffolk to the south, Cambridgeshire and Essex to the west.

A.5 Institution ID

Please contact DCC if your institution is not listed, or you wish to modify the text

ID	Institution				
001	Spanish National Cancer Research Centre (CNIO, Madrid)				
002	Barcelona Supercomputing Center (BSC-CNS, Madrid)				
003	EMBL-EBI (Hinxton)				

Table A.6.1: Value Codes for Platform or Validation Platform

A.6 Value Codes for Platform or Validation Platform

Please contact the DCC if your platform/technology is not listed here.

Key	Platform or Validation Platform
1	PCR
2	qPCR
3	capillary sequencing
4	SOLiD sequencing
5	Illumina GA sequencing
6	454 sequencing
7	Helicos sequencing
8	Affymetrix Genome-Wide Human SNP Array 6.0
9	Affymetrix Genome-Wide Human SNP Array 5.0
10	Affymetrix Mapping 100K Array Set
11	Affymetrix Mapping 500K Array Set
12	Affymetrix Mapping 10K 2.0 Array Set
13	Affymetrix EMET Plus Premier Pack
14	Agilent Whole Human Genome Oligo Microarray Kit
15	Agilent Human Genome 244A
16	Agilent Human Genome 105A
17	Agilent Human CNV Association 2x105K
18	Agilent Human Genome 44K
19	Agilent Human CGH 1x1M
20	Agilent Human CGH 2x400K
21	Agilent Human CGH 4x180K
22	Agilent Human CGH 8x60K
23	Agilent Human CNV 2x400K
24	Agilent Human miRNA Microarray Kit (v2)

Agilent Human CpG Island Microarray Kit

Agilent Human SpliceArray Illumina human1m-duo

Agilent Human Promoter ChIP-on-chip Microarray

Continued on next page

26

27

Set

Table A.6.1 – continued from previous page

v3.0
2.0D
nome
nome
e Set
07-1
07-2
07-3
oaary
page

Table A.6.1 – concluded from previous page

Key	Platform or Validation Platform
67	Ion Torrent PGM
68	Illumina GoldenGate Methylation Cancer Panel I
69	Illumina Infinium HumanMethylation450
70	Agilent 8 x 15K Human miRNA-specific microarray
71	M.D. Anderson Reverse Phase Protein Array Core
7 2	Microsatellite Instability Analysis
7 3	Agilent 244K Custom Gene Expression G4502A-07
74	Illumina HumanCNV370-Duo v1.0 BeadChip
75	Illumina HumanOmniExpress BeadChip

A.7 Chromosome Names for Reference Genome GRCh37

Table A.7.1: Chromosome Names for Reference Genome GRCh37

Key	Chromosome Name
1	1
2	2
3	3
4	4
5	5
6	6
7	7
8	8
9	9
10	10
11	11
12	12
13	13
14	14
15	15
16	16
17	17
18	18
19	19
20	20
21	21
22	22
23	X
24	Υ
25	MT
26	c5_H2
	Continued on next page

Table A.7.1 – continued from previous page

	Table A.7.1 – continued from previous page
Key	Chromosome Name
27	c6_COX
28	c6_QBL
29	NT_113870
30	NT_113871
31	NT_113872
32	NT_113874
33	 NT_113878
34	NT_113880
35	NT_113881
36	NT_113884
37	NT_113885
38	NT_113886
39	NT_113888
40	NT_113889
41	NT_113890
42	NT_113898
43	NT_113899
44	NT_113901
45	NT_113902
46	NT_113903
47	NT_113906
48	NT_113908
49	NT_113909
50	NT_113910
51	NT_113911
52	NT_113912
53	NT_113915
54	NT_113916
55	NT_113917
56	NT_113923
57	NT_113924
58	NT_113925
59	NT_113926
60	NT_113927
61	NT_113929
62	NT_113930
63	NT_113931
64	NT_113932
65	NT_113933
66	NT_113934
67	NT_113935
68	NT_113936
	Continued on next page

Table A.7.1 – continued from previous page

	Table A.7.1 – continued from previous page
Key	Chromosome Name
69	NT_113937
70	NT_113939
71	NT_113943
7 2	NT_113944
73	NT_113946
74	NT_113949
75	NT_113951
76	NT_113953
77	NT_113954
78	NT_113956
79	NT_113957
80	NT_113958
81	NT_113960
82	NT_113961
83	NT_113962
84	NT_113963
85	NT_113964
86	NT_113965
87	NT_113966
88	HSCHR17_1
89	HSCHR17_RANDOM_CTG2
90	HSCHR17_RANDOM_CTG3
91	HSCHR19_RANDOM_CTG2
92	HSCHR1_RANDOM_CTG12
93	HSCHR1_RANDOM_CTG5
94	HSCHR4_RANDOM_CTG2
95	HSCHR4_RANDOM_CTG3
96	HSCHR6_MHC_APD
97	HSCHR6_MHC_COX
98	HSCHR6_MHC_DBB
99	HSCHR6_MHC_MANN
100	HSCHR6_MHC_MCF
101	HSCHR6_MHC_QBL
102	HSCHR6_MHC_SSTO
103	HSCHR7_RANDOM_CTG1
104	HSCHR8_RANDOM_CTG1
105	HSCHR8_RANDOM_CTG4
106	HSCHR9_RANDOM_CTG2
107	HSCHR9_RANDOM_CTG4
108	HSCHR9_RANDOM_CTG5
109	HSCHRUN_RANDOM_CTG1
110	HSCHRUN_RANDOM_CTG10
	Continued on next page

Table A.7.1 – concluded from previous page

1/	Table 71.7.1 – concluded from previous page
Key	Chromosome Name
111	HSCHRUN_RANDOM_CTG11
112	HSCHRUN_RANDOM_CTG13
113	HSCHRUN_RANDOM_CTG14
114	HSCHRUN_RANDOM_CTG15
115	HSCHRUN_RANDOM_CTG16
116	HSCHRUN_RANDOM_CTG17
117	HSCHRUN_RANDOM_CTG2
118	HSCHRUN_RANDOM_CTG20
119	HSCHRUN_RANDOM_CTG21
120	HSCHRUN_RANDOM_CTG22
121	HSCHRUN_RANDOM_CTG23
122	HSCHRUN_RANDOM_CTG26
123	HSCHRUN_RANDOM_CTG29
124	HSCHRUN_RANDOM_CTG3
125	HSCHRUN_RANDOM_CTG30
126	HSCHRUN_RANDOM_CTG31
127	HSCHRUN_RANDOM_CTG32
128	HSCHRUN_RANDOM_CTG33
129	HSCHRUN_RANDOM_CTG34
130	HSCHRUN_RANDOM_CTG35
131	HSCHRUN_RANDOM_CTG36
132	HSCHRUN_RANDOM_CTG4
133	HSCHRUN_RANDOM_CTG40
134	HSCHRUN_RANDOM_CTG5
135	HSCHRUN_RANDOM_CTG6
136	HSCHRUN_RANDOM_CTG9
137	HSCHR4_1

A.8 Value Codes for Reference Genome Assembly Version

Table A.8.1: Value Codes for Reference Genome Assembly Version

Key	Reference Genome Assembly Version
1	GRCh37
2	NCBI36
3	GRCh37.p1
4	GRCh37.p2
5	GRCh37.p3
6	GRCh37.p4
7	GRCh37.p5

A.9 Value Codes for Raw Data Repository

Table A.9.1: Value Codes for Raw Data Repository

Key Raw Data Repository

Continued on next page

Table A.9.1 – concluded from previous page

Key	Raw Data Repository
1	EGA
2	dbSNP
3	TCGA
4	CGHub
5	GEO